


SMT5510 Scope container - Container with silicone dividers and lid

Long Life Re-usable Polypropylene Container for Scopes and Long Instruments for use in Sterile Services and Decontamination Units

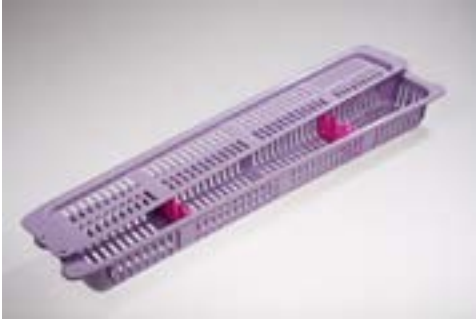

The Warwick Range of polypropylene products from Warwick SASCo was originally designed to offer an alternative to stainless steel products traditionally used in Hospitals.

This Range is now widely used in Sterile Services and Decontamination Units, Operating Rooms and Clinics due to its low cost, effective product design, and repeat performance results in hospital autoclaves.

Easily recognised by its distinct blue colour, the Warwick Range includes bowls, gallipots, kidney dishes, instrument trays, instrument protection, instrument boxes, traceability tags, measuring jugs, scope and instrument soaking containers.



Product Benefits

- Mesh base, sides and lid allow maximum water penetration and steam circulation.
- Moveable silicone grips hold delicate scopes and instruments in place - reducing damage to instruments during cleaning, sterilization and transportation.
- Lightweight, strong and durable - 50% lighter than stainless steel.
- Improved removal of blood residue compared to stainless steel.
- Holes for traceability and security tags.
- Rapid water draining gives improved drying results due to raised floor and internal curves.
- Stackable design for easy storage.
- Smooth corners and edges prevent damage to wrapping papers.
- Maximise Sterile Services production as tray saves space during processing.
- Can be sterilized by Autoclave, Ethylene Oxide or Hydrogen Plasma.

Product Specification

- Mesh Tank with Lid and Silicone Grips
- External dimensions: 660mm x 110mm x 66mm
- Internal dimensions: 600mm x 85mm x 55mm
- Pack Size 1
- The following standards are applicable:
 - **ISO 9001:2008**
Quality Management System
 - **ISO 13485:2003**
Medical Devices - Quality Management System
 - **ISO 17664:2004**
Sterilization of Medical Devices - Guidelines for the processing of re-sterilizable medical devices
- Manufactured from high quality polypropylene
- Autoclavable up to 135 ° C

Guide to Using Warwick Products from Warwick SASCo

Cleaning

- Product should be washed before first use.
- Warwick Products are designed for a lengthy service and guaranteed for 1000 washes and sterilizations. Scouring and abrasive cleaners will reduce the performance of products.
- Product should be processed complying with requirements of ISO17664.

Storage

- This product has a shelf life of 12 months from date of manufacture providing the recommended storage conditions are adhered to.
- Product should be stored in its primary packaging, off the floor on a pallet or shelf, and in a clean dry ventilated room at a temperature of 4 ° C - 25 ° C.